

User Manual



Smart Mini Massage Gun GS07

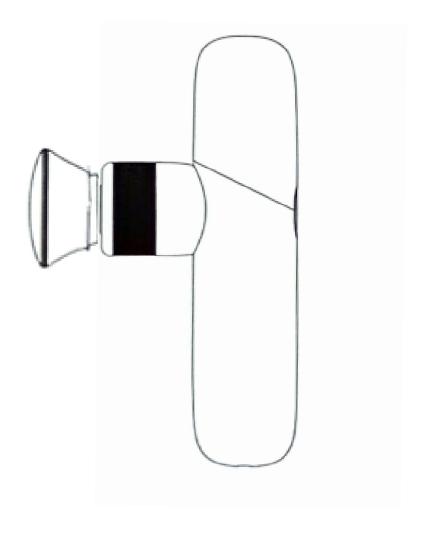

Product Introduction

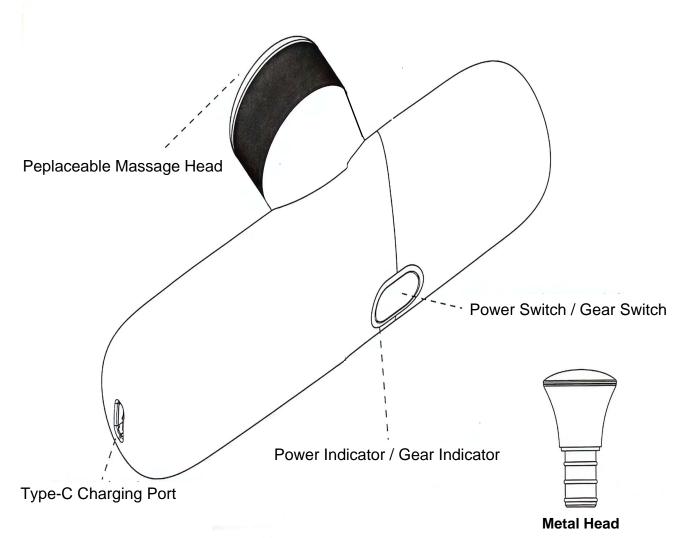

Massage Head

Suitable for massaging arm, waist, buttocks, thighs, calves and other muscle groups.



U-Shaped Head

and achilles tendon.



Cone Head

Suitable for neck, spine Suitable for impacting deep tissues such as meridians, palms,oles, etc.



Round Head


Suitable for arm, waist, buttocks, thighs, calves and other muscle groups.




Flat Head

Suitable for relaxation and shaping of all parts.


Installing / Replacing Massage Heads

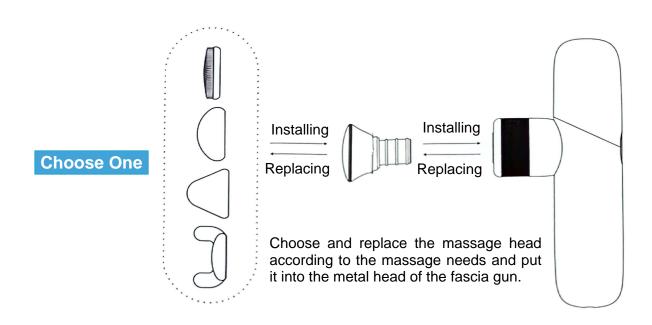

Installation Effect

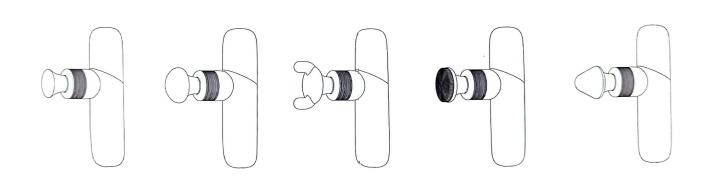


On/Off

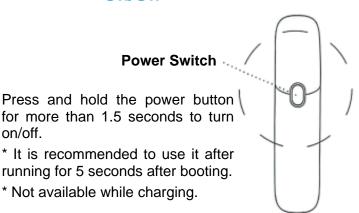

Gear Switching

After powering on, it enters the smart mode, and the color of the indicator light is displayed according to the user's different usage strengths: white, blue, green, and red.


Short press 2 times after powering on to enter normal mode; 1,2, 3, and 4 gears can be switched in a cycle.

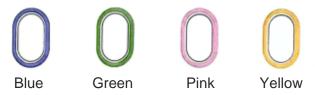

- * Because it is a metal body, it will heat up after prolonged use, which is a normal phenomenon.
- * Excessive pressing force will cause the machine to stop, all gear lights will flash, short press the power button to start running.
- * To protect the motor, the machine will enter the protection state (automatic shutdown) after 10 minutes of use. If you need to continue to use it, you can restart it.

Battery Level Indication


In the power-off state, short press the power button to display the battery level, and the color light will turn off for 1.5 seconds.

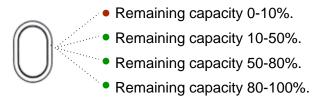


Type-C Charging Port Type-C Charging Port SV/2A Power Indicator While charging • Always bright-charging. While charging • Always bright-end of charge.

*This massage gun does not support the fast charging function of mobile phone.

Packaging List

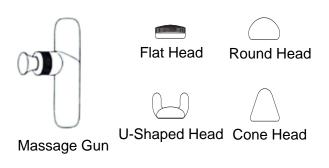

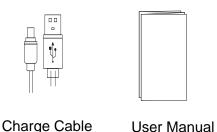

Specifications:

Shock stroke: 7.5mm Thrust: >17 lbs

Battery capacity: 1000mAh Noise: <45dB (1 gear)

Charging time: <3 hours Product size: 5.7x2.6x1.4 IN
Battery time: >4 hours Machine net weight: 0.6 lbs

Charging port: USB-C Gross weight: 0.8 lbs

Input voltage: DC 5V1.2A

Store and Save

Please store this product out of the reach of children and pets.

- Do not place this product near a heat source (stove or heater, etc.) and in a car on a hot day.
- Please store in a dry environment, do not overcharge or over-discharge the battery, otherwise it will cause damage to the battery.
- Do not use this product in an environment where the temperature is too high or too low.
- The polymer battery used in this product is an important resource that can be recycled and regenerated. When the product is discarded, be sure to disassemble the battery and hand it over to a professional organization with battery recycling capabilities for disposal.

Disclaimer

Read the instructions carefully before using the massager. The use of this equipment is subject to following its instructions closely. ABEGET is not liable for any damage or loss caused by inappropriate device use.

Do not use the massager without the approval of your doctor if you are: pregnant, diabetic, have a pacemaker, recovering from surgery, suffer from migraines, neuropathy disorders, disc protrusion, or spinal inflammation.

If you feel pain or discomfort while using the massager, please stop using it immediately and consult your doctor or physiotherapist.